

The Carl-Oskar is designed for users that find it difficult to maneuver a regular rollator or require extra support. The forearm supports are adjustable in height and angle so that the user finds the optimal positioning for comfort and stability.

#### **Outstanding Stability**

The Carl-Oskar RA is the ideal rollator for users recovering from stroke. It features separate forearm supports and large ball-bearing mounted wheels that ensure walking is easy, even on uneven or rough surfaces. Its reactive frame is durable, even for the more active users. The Carl-Oskar RA is available with a seat as an option.

#### **Unique Handles**

The Carl-Oskar RA features unique front-facing, angle adjustable handles that allow users to find the optimal handle position for comfort and stability. Adjustment is easy with integrated measurements right in the handle tube. The ergonomic design of the

grips ensures even distribution of weight. The Carl-Oskar RA can also be equipped with a one-hand brake to allow a user to control hand functions on either the left or right hand side.

#### **Adjustable Forearm Supports**

The Carl-Oskar RA features very unique forearm supports that are adjustable in both height and direction. This ensures that the forearm supports can be positioned in the best possible way for user comfort and stability.

#### **HUMAN CARE SWEDEN (HQ)**

Årstaängsvägen 21B  
117 43 Stockholm  
Phone:- +46 8 665 35 00  
Fax:- +46 8 665 35 10  
info.se@humancaregroup.com

#### **HUMAN CARE CANADA**

10-155 Colonnade Road  
Ottawa, ON K2E 7K1  
Phone: +613.723.6734  
Fax: +613.723.1058  
info.ca@humancaregroup.com

#### **HUMAN CARE NEDERLAND**

Elspeterweg 124  
8076 PA Vierhouten  
Phone: +31 577 412 171  
Fax: +31 577 412 170  
info.nl@humancaregroup.com

#### **HUMAN CARE UNITED STATES**


8006 Cameron Road, Suite K  
Austin, TX 78754  
Phone: +512.476.7199  
Fax: +512.476.7190  
info.us@humancaregroup.com

## Accessories/Features



### Tray

Practical for resting lightweight items. Maximum load 2kg/4.4lb. Part no. 10601

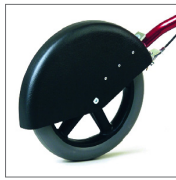

### Splash Guard

Protects legs and clothes. Part no. 10608



### Slow-Down Brake

Prevents the rollator from moving too fast and enables a controlled walking speed for the user. Part no.: 10624



### One-Hand Brake

Allows user to control hand functions with one hand (installed left or right). Special order item. Part no. 10653




### Cane holder

Simplifies transport of a cane or umbrella. Part no. 10603



### Curb Climber

Helps raise rollator onto a curb. Part no. 10612

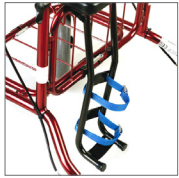

### Oxygen Holder

Holds oxygen tube firmly in place. Part no. 10615



### Basket

Included with rollator in Canada. Ideal for shopping Part no. 10601

## Technical Specifications

	Carl-Oskar RA Low 67	Carl-Oskar RA Low 72	Carl-Oskar RA Nor 67	Carl-Oskar RA Nor 72
Part no.	12315	12415	12325	12425
Part no. with seat	12315-1	12415-1	12325-1	12425-1
Handle Height (Min - Max)	95cm - 113cm/ 37.4" - 44.5"	95cm - 113cm/ 37.4" - 44.5"	107cm - 126cm/ 42.1" - 49.6"	107cm - 126cm/ 42.1" - 49.6"
Seat Height	64cm/25.2"	64cm/25.2"	64cm/25.2"	64cm/25.2"
Overall Width	67cm/26.4"	72cm/28.3"	67cm/26.4"	72cm/28.3"
Width between Handles	47cm/18.5"	52cm/20.5"	47cm/18.5"	52cm/20.5"
Distance between Rear Wheels	52.1cm/20.5"	52.1cm/20.5"	52.1cm/20.5"	52.1cm/20.5"
Wheel Size	30cm x 5cm/ 11.8" x 1.8"	30cm x 5cm/ 11.8" x 1.8"	30cm x 5cm/ 11.8" x 1.8"	30cm x 5cm/ 11.8" x 1.8"
Product Weight	13kg/28.6lb	13.1kg/28.9lb	13.2kg/29.1lb	13.3kg/29.3lb
Maximum User Weight	125kg/275.6lb	125kg/275.6lb	125kg/275.6lb	125kg/275.6lb

Note: Not all accessories are available for sale in all countries.