

Comfort Sling Soft & Comfort Sling High Soft



Sit on Sling

Undivided Legs

Comfortable



Comfort Sling Soft High
Art. Nr. 25069



Comfort Sling Soft
Art. Nr. 25064

Our Comfort Slings are designed to give the patient comfort and a sense of safety. The undivided legs evenly distribute the patient's weight in order to reduce pressure, making this a suitable option for patients with pain. These slings come in two models, one with extra head support.

Unique Fabric

Our Comfort Slings are made of Soft Air Mesh, a material exclusive to Human Care. This soft material promotes less repositioning making it the optimal choice for patients with pain. In addition, the material does not promote skin breakdown and can therefore be left behind the patient for short periods of time.

Functional Design

Our Comfort Slings have no reinforcements or handles and feature undivided legs that evenly distribute the patient's weight. With even distribution, pressure is reduced and patient comfort increased. Our Comfort Slings are available in a medium back and high back

model that offers extra head support. Both slings are considered "sit-on" slings and are applied when the patient is lying down.

Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.



Walking Aids



Healthcare Beds



Lifting Solutions



Hygiene



Convertible Chairs

Accessories



Extension Loops
Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no.:21180-*



Leg Padding and Supports
Ideal for users with sensitive skin or pain. Part no.: 21146-S, M, L Reinforced supports. One size. Part no. 25102



Upholstered Cover
Protects user from hanger bar impact. Not for use in a wet environment. Part no. 50866



Hanger Bars for 2-, 3- or 4-point Suspension
Part no. 57551, 55764 (2- or 3-Point - Quick Release), 55619 (4-Point - Attach to leg hook)


Lifting Situation




Use to/from a bed
Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



Use to/from a chair
Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.

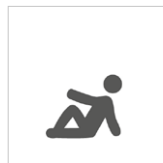


Use to/from the floor
Suitable for transfers from the floor.

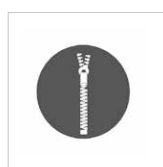


Use for dressing/undressing
Suitable for dressing/undressing applications.

Care and Features

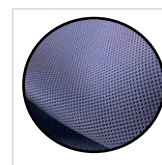

Soft Air Mesh
Lightweight, cool and breathable material. Suitable for sensitive skin.

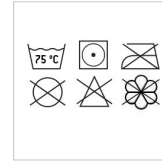

Care
Care instructions located on sling label. Use a hospital grade laundry bag while laundering

Technical specifications

Part Number	Size	Colour	Model #25064			Model #25069				
			A	B	C	A	B	C		
Comfort Sling Soft	Comfort Sing High Soft									
25064-JM	25069-JM	Junior Medium		76cm / 30"	73cm / 28.7"	90cm / 35.4"	76cm / 30"	73cm / 28.7"	111cm / 43.7"	
25064-S	25069-S	Small		76cm / 30"	80cm / 31.5"	100cm / 39.3"	76cm / 30"	80cm / 31.5"	122cm / 48"	
25064-M	25069-M	Medium		86cm / 33.9"	90cm / 35.4"	110cm / 43.3"	86cm / 33.9"	90cm / 35.4"	134cm / 52.8"	
25064-L	25069-L	Large		96cm / 37.8"	105cm / 41.3"	116cm / 45.7"	96cm / 37.8"	105cm / 41.3"	139cm / 54.7"	
25064-XL	25069-XL	Extra Large		106cm / 41.7"	118cm / 46.5"	120cm / 47.2"	106cm / 41.7"	118cm / 46.5"	145cm / 57.1"	
Maximum User Weight	300kg / 660 lb									

* Available in assorted lengths: 10cm/3.9" (21180-10), 15cm/5.9" (21180-15), 20cm/7.8" (21180-20), 25cm/9.8" (21180-25), 35cm/13.7" (21180-35)



We improve freedom of movement and quality of care for people.

www.humancaregroup.com

