



# Toileting Sling & Toileting Sling High



Padded Arms  
Reinforced Back  
Divided Legs

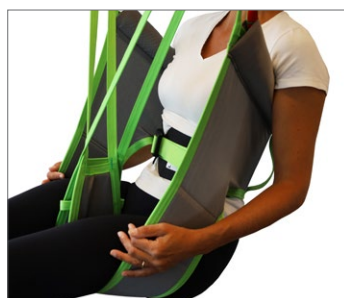

Art. nr. 25080B, Toileting Sling w. belt

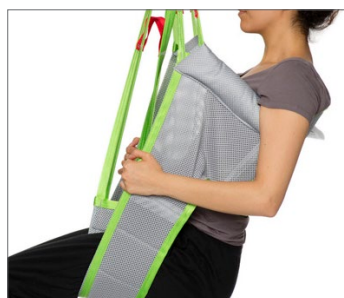

Art. nr. 25085, Toileting Sling Net



Art. nr. 25090, Toileting Sling High Net

Our Toileting Slings feature a wide, divided leg partition that enables fast and simple application for dressing and undressing during toileting or bathing. They are suitable for patients with poor trunk stability and feature padding under the arms for extra comfort and are available with extra head support or without.

### Three Options

These Toileting slings are available in two low back models and one high back model with adjustable neck support for increased head stability. The wide partition of the leg parts makes the sling application fast and simple and enables dressing and undressing during toileting or bathing. All models can be used with or without an optional waist belt that is available in an assortment of sizes with a Velcro or buckle closure.

### Optimal Positioning

The low back models feature an upright seated position that makes transferring to and from a toilet much easier. The high back model features a reclined (backward tilted)

position that prevents the user from sliding out of the sling and is especially beneficial for those hypertonic users that have a tendency to overstress themselves (spastic performance) or those hypotonic users with lesser body control.

### Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.



Walking Aids



Healthcare Beds



Lifting Solutions



Hygiene



Convertible chairs

## Accessories



### Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no. 21180.\*



### Waist Belts

Velcro or buckle enclosure. Available in Junior and Adult sizing. Velcro Part no. 25145-\*\* Buckle Part no. 25146-\*\*\*



### Upholstered Cover

Protects user from hanger bar impact. Not for use in a wet environment. Part no. 50866



### Hanger Bars for 2-, 3- or 4-point Suspension

Part no. 57551, 55764 [2- or 3-Point - Quick Release], 55619 [4-Point - Attach to leg hook]

## Lifting Situations



### Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



### Use to/from a chair

Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



### Use for showering


Only Polyester Net fabric can be used for this application.



### Use for toileting


Only Polyester Net fabric can be used for this application.

## Care and Features

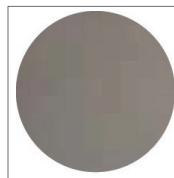


### Polyester

Durable, low-friction, easy to clean, apply and remove.

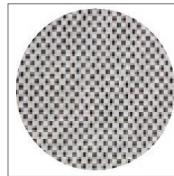

### Polyest Net

Soft, thin material that breathes and dries quickly. Ideal for bathing or showering.

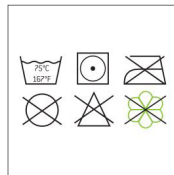

### Care

Care instructions located on sling label. Use a hospital grade laundry bag while laundering.



## Technical Specifications

Article number		Size	Colour	Model #25080/25085			Model #25090			Modell #25080/25085	Modell #25090
				A	B	C	A	B	C		
Toileting Sling Net / Toileting Sling High Net	Toileting Sling Polyester	Toileting sling with belt									
-	25080-JS	25080B-JS	Junior Small	JS	28cm / 11"	62cm / 24.4"	22cm / 8.7"	-	-	-	
25085-JM / 25090-JM	25080-JM	25080B-JM / 25090B-JM	Junior Medium	JM	28cm / 11"	67cm / 26.3"	26cm / 10.2"	28cm / 11"	73cm / 28.7"	70cm / 27.6"	
25085-S / 25090-S	25080-S	25080B-S / 25090B-S	Small	S	32cm / 12.6"	74cm / 29.1"	30cm / 11.8"	32cm / 12.6"	82cm / 32.3"	74cm / 29.1"	
25085-M / 25090-M	25080-M	25080B-M / 25090B-M	Medium	M	48cm / 18.9"	90cm / 35.4"	33cm / 13"	48cm / 18.9"	92cm / 36.2"	83cm / 32.7"	
25085-L / 25090-L	25080-L	25080B-L / 25090B-L	Large	L	49cm / 19.3"	97cm / 38.2"	33cm / 13"	49cm / 19.3"	99cm / 39"	90cm / 35.4"	
25085-XL / 25090-XL	25080-XL	25080B-XL / 25090B-XL	Extra Large	XL	56cm / 22"	100cm / 39.4"	38cm / 15"	56cm / 22"	100cm / 39.4"	100cm / 39.4"	
Högsta användarvikt: 300 kg / 660 lb											

\* Available in assorted lengths: 10cm/3.9" [21180-10], 15cm/5.9" [21180-15], 20cm/7.8" [21180-20], 25cm/9.8" [21180-25], 35cm/13.7" [21180-35]

\*\* Junior Sizes: 40cm/15.7" [25145-40B], 50cm/19.7" [25145-50B], 60cm/23.6" [25145-60B]; Adult Sizes: 60cm/23.6"/[25145-60], 70cm/27.6" [25145-70], 80cm/31.5" [25145-80], 90cm/35.4" [25145-90], 110cm/43.3" [25145-110], 120cm/47.3" [25145-120]

\*\*\* Junior Sizes: 40cm/15.7" [25146-40B], 50cm/19.7" [25146-50B], 60cm/23.6" [25146-60B]; Adult Sizes: 60cm/23.6"/[25146-60], 70cm/27.6" [25146-70], 80cm/31.5" [25146-80], 90cm/35.4" [25146-90], 100cm/39.4" [25146-100], 110cm/43.3" [25146-110], 120cm/47.3" [25146-120]



We improve freedom of movement and quality of care

www.humancaregroup.com

**HUMAN CARE SWEDEN (HQ)**  
Årstaängsvägen 21B  
117 43 Stockholm  
+46 8 665 35 00  
info@humancaregroup.com

**HUMAN CARE CANADA**  
10-155 Colonnade Road  
ON K2E 7K1 Ottawa  
+1 613 723 6734  
info.ca@humancaregroup.com

**HUMAN CARE THE NETHERLANDS**  
Elsperweg 124  
8076 PA Vierhouten  
+31 577 41 21 71  
info.nl@humancaregroup.com

**HUMAN CARE UNITED STATES**  
8006 Cameron rd, Suite K  
TX 78754 Austin  
+1 5 124 76 7199  
info.us@humancaregroup.com