


Positioning Slings



Ventilated
Bariatric Model
SPU Option



Our positioning slings are available in reusable or single patient use models and are designed to assist a caregiver to safely and securely perform safe patient handling during sudden or intended positioning/repositioning, turning and lateral transfer with a patient in the supine position.

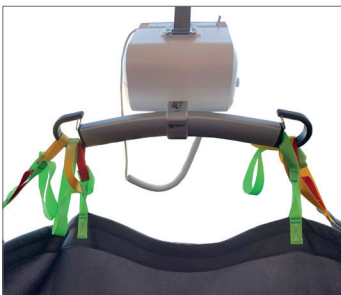


Ventilated Material

Our reusable ventilated slings are made of a breathable, plyable, material that has moisture-wicking properties. The unique material facilitates air flow making these slings ideal for critical care units as they can be placed behind the patient like a sheet for sudden and often unpredictable repositioning requirements. These slings easily fit most hospital beds and the material does not promote skin breakdown allowing the sling to remain behind the patient for short periods of time.

the added benefit of hygiene control. Made of a durable Polyester/Cotton blend, or Soft Air Mesh and Derma Therapy material - these slings can be used by the patient for multiple lifts before disposal.

Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature eight reinforced straps (6 body/2 head) and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. The straps are positioned to ensure optimal comfort and safety to the patient.

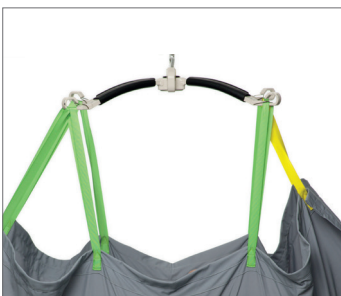

Single Patient Use

Our single patient use slings offer the convenient option of disposal after use with

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone: +46 8 665 35 00
Fax: +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com


HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES


8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories

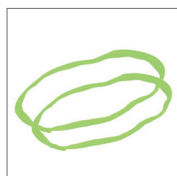


Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no.:21180.*

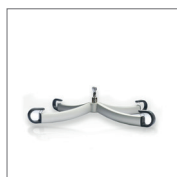

Hanger Bar for 4-point Suspension

92938 (4-Point for HeliQ or HeliQ PT)

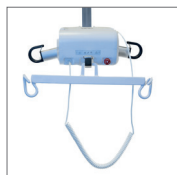

Hanger Bar for 4-point Suspension

55619 (4-Point - Attach to leg hook on Altair/Roomer S)

Lifting Situations



Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



Use to/from a chair


Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



Single Patient Use Model


Dispose after use. Intended for short term use by a single patient.

Care and Features

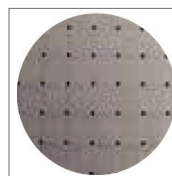


Ventilated Polyester

Maximum airflow and breathability. Can be left behind patient like a sheet.

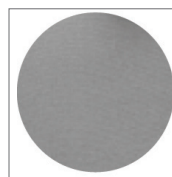


Polyester/Cotton

Single patient use.

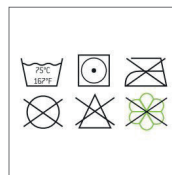


Care - Reusable Slings

Care instructions located on sling label. Use a hospital grade laundry bag while laundering.

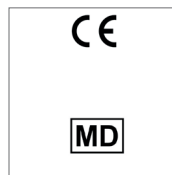

Technical Specifications

Part Number	Weight Capacity	A	B
Reusable Polyester, Ventilated Soft Air Mesh	SWL kg/lb	Width cm/inch	Length cm/inch
NSB-410-ST-L	454kg/1000lb	1320mm / 52"	208cm/81.9"
NSB-410-ST-XL	454kg/1000lb	147cm/57.9"	208cm/81.9"
NSB-451-S	318 kg/700lb	1050 mm/41.3"	1680 mm/66"
Single Patient Use Polyester/Cotton Soft Air Mesh and DermaTherapy			
NSB-451-U	318kg/700lb	105cm /41.3"	168cm/ 66.1"
NSB-451-A			
NSB-451-DT	385 kg/850lb	1090 mm /43"	2180 mm/85,8"

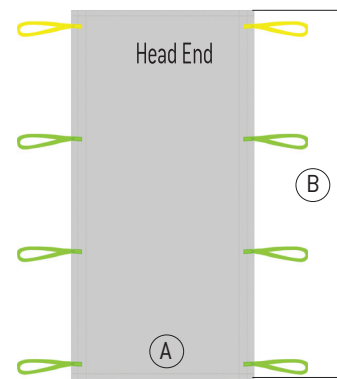

* Available in assorted lengths: 10cm/3.9" (21180-10), 15cm/5.9" (21180-15), 20cm/7.8" (21180-20), 25cm/9.8" (21180-25), 35cm/13.7" (21180-35)

** Only available in the United States