



Safe E-clips

Ergonomic support

Comfortable fit



Support during standing up

When helping a client stand up, care givers often struggle to keep things safe and ergonomic for both people involved. Our Back Belt E-Clip provides just the support they may need to properly use our seBastian II sit-to-stand lift.



The **Back Belt** encourages natural movement, with the padded back providing comfortable back support while getting up. The Back Belt also ensures the client feels safe and secure and prevents any risk of falling backwards. The integrated waist belt with buckle fastening keeps the sling in place during the entire transfer. The non-slip surface on the inside of the back prevents the sling from sliding up while the client gets up, allowing both the client and the care giver to focus entirely on standing up. For clients with limited core stability or who have other reasons that may make getting up difficult, the Back Belt E-clip is a welcome addition when using a sit-to-stand lift.

The **universal E-clip** is a proven system for fastening a sling to a lift, thereby significantly reducing the risk of error. The clip is placed around the snap hooks of the hanger/spreader bar and closes with an audible 'click'. It is a signal to the care giver that the clip is securely attached and cannot unexpectedly disconnect.

Our slings are easy to fit, clean and maintain and are colour-coded to help you find the right size quickly and easily. Human Care slings are produced to the highest quality standards. The sling is surrounded by a strong green harness for excellent strength and durability.



Mobility Aids



Healthcare Beds



Lifting Solutions



Hygiene



Convertible Chairs

Use



From and into bed

Suitable for lifting clients from bed and helping them sit down on bed.




From and to a chair

Suitable for positioning in a chair and for transfers in and out of a chair.



Standing exercises

Suitable for standing and rehabilitation exercises.

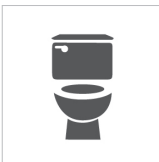

Using the toilet

Provides extra support when using the toilet.

Features




Polyester

Durable material with a low coefficient of friction, making it easy to apply and clean.



E-clip integrated lock

Can be used one-handed. Universally compatible with all lifts equipped with mushroom fasteners Ø 6.8-8.2 mm.

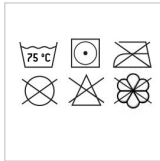


Cleaning and maintenance instructions

For cleaning and care instructions, see the label of the sling. Wash in an approved laundry bag.



Technical specifications

Product name	Article number	Size	Colour	A	B	C
Back Belt E-clip S	25130C-S	Small		69 cm / 27.2"	27 cm / 10.6"	33cm / 13"
Back Belt E-clip M	25130C-M	Medium		83 cm / 32.7"	30 cm / 11.8"	39 cm / 15.3"
Back Belt E-clip L	25130C-L	Large		92 cm / 36.2"	34 cm / 13.4"	44 cm / 17.3"
Back Belt E-clip XL	25130C-XL	Extra Large		103 cm / 40.5"	39 cm / 15.3"	50 cm / 19.7"
	Maximum user weight	250kg / 551lb				

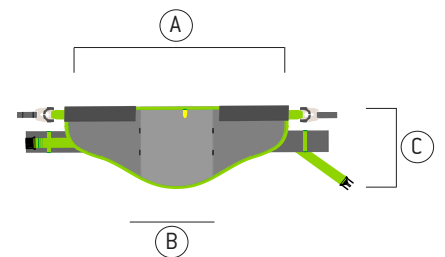

We improve freedom of movement and quality of care for people.

www.humancaregroup.com



HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
+46 8 665 35 00
info@humancaregroup.com

HUMAN CARE GERMANY

Industriestraße 13, Gebäude 1/9
63755 Alzenau
+49 (0) 6023 994 936 8
infode@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
ON K2E 7K1 Ottawa
+1 613 723 6734
info.ca@humancaregroup.com

HUMAN CARE THE NETHERLANDS

Els peterweg 124
8076 PA Vierhouten
+31 577 41 21 71
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron rd, Suite K
TX 78754 Austin
+1 5 124 76 7199
info.us@humancaregroup.com