




Pro Skin Positioning Sling

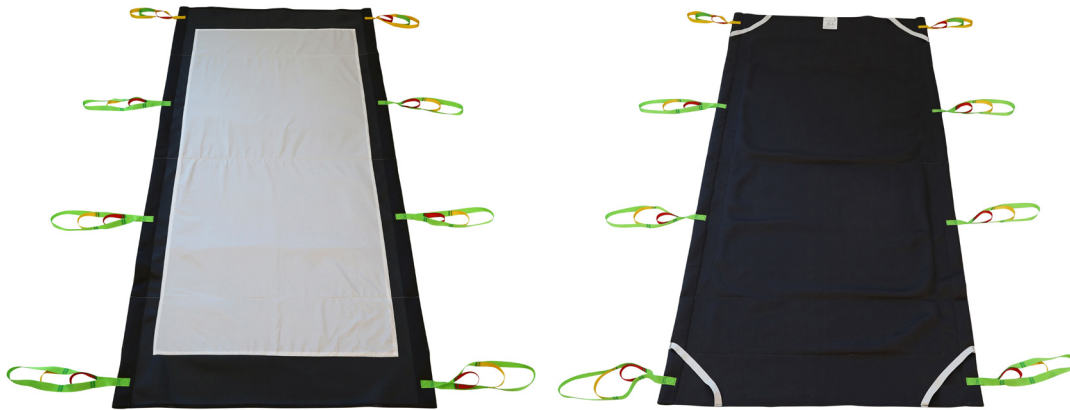

Pro-skin
Gentle
Positioning

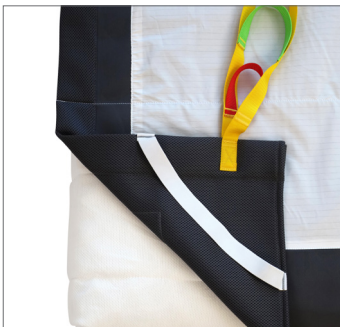

The unique design of the Pro Skin Positioning Sling incorporates DermaTherapy material which decreases the risk of pressure ulcers in clients. The sling has elastic bands on each corner in order to better fit onto the clients' bed. The eight loops are color coordinated with yellow loops by the the head end of the bed for guidance.

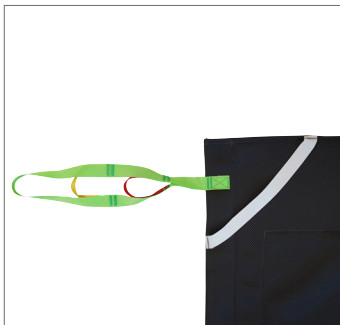

Human Care positioning slings

Our positioning slings are available in reusable or single patient use models. They are designed to assist a caregiver to safely and to securely perform safe patient handling during intended or sudden positioning/repositioning, turning and lateral transfer with a patient in the supine position.

High quality material

The Soft Air Mesh material allows moisture and heat to transfer efficiently from the patient which allows the sling to be left behind the patient for a period of time after the transfer is complete. The DermaTherapy® material which is incorporated into this positioning sling creates a cleaner, drier and hypoallergenic experience. The DermaTherapy material efficiently decreases the risk for pressure ulcers.



Used together with a 2-point hanger bar carers are enabled to reposition a client in bed; over to the side or up and down. Used together with a 4-point hanger bar the sling can also assist transfers during situations like 'bed to bed', from 'floor to bed' and from 'bed into a seated position in chair'.

Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. The Pro Skin Positioning Sling features eight reinforced loops (6 body/2 head) and high strength nylon webbing around the sling to make them strong enough to resist even the most active



Walking Aids



Healthcare Beds



Lifting Solutions



Hygiene



Convertible chairs

Accessories



Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no.:21180-*



Hanger Bar for 4-point Suspension

4-Point - for HeliQ or HeliQ PT
Article number: 92938



Hanger Bar for 4-point Suspension

4-Point - Attach to leg hook on Altair/Roomer S
Article number: 55619

Lifting Situations



Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



Use to/from a chair

Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



Single Patient Use Model


Dispose after use.
Intended for short term use by a single patient.



Use to/from the floor

Suitable for transfers from the floor.

Care and Features

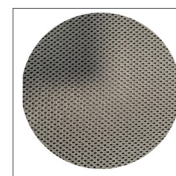

Soft Air Mesh

Lightweight, cool and breathable material. Suitable for sensitive skin.



Technical Specifications

Part Number	Weight Capacity	A	B	
NSB-451-DT	385 kg / 850 lbs	109 cm / 42,9"	218 cm / 85,8"	<p>Head end</p> <p>Yellow loops by head end</p>

www.humancaregroup.com

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone:- +46 8 665 35 00
Fax:- +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com